



HAL
open science

A multiparameter optimization in middle-down analysis of monoclonal antibodies by LC–MS/MS

Jonathan Dhenin, Mathieu Dupré, Karen Druart, Alain Krick, Christine Mauriac, Julia Chamot-Rooke

► **To cite this version:**

Jonathan Dhenin, Mathieu Dupré, Karen Druart, Alain Krick, Christine Mauriac, et al.. A multiparameter optimization in middle-down analysis of monoclonal antibodies by LC–MS/MS. *Journal of Mass Spectrometry*, 2023, 58 (3), pp.e4909. 10.1002/jms.4909 . pasteur-04102833

HAL Id: pasteur-04102833

<https://pasteur.hal.science/pasteur-04102833>

Submitted on 22 May 2023

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons Attribution - NonCommercial - NoDerivatives 4.0 International License

A multiparameter optimization in middle-down analysis of monoclonal antibodies by LC-MS/MS

Jonathan Dhenin^{1,2,3}  | Mathieu Dupré¹ | Karen Druart¹ | Alain Krick³ |
Christine Mauriac³ | Julia Chamot-Rooke¹

¹Institut Pasteur, Université Paris Cité, CNRS UAR2024, Mass Spectrometry for Biology, Paris, 75015, France

²Université Paris Cité, Sorbonne Paris Cité, Paris, France

³DMPK, Sanofi, Chilly-Mazarin, 91385, France

Correspondence

Julia Chamot-Rooke, Institut Pasteur, Université Paris Cité, CNRS UAR2024, Mass Spectrometry for Biology, F-75015 Paris, France.

Email: julia.chamot-rooke@pasteur.fr

Jonathan Dhenin, DMPK, Sanofi, Chilly-Mazarin, France.

Email: jonathan.dhenin@sanofi.com

Funding information

Association Nationale de la Recherche et de la Technologie, Grant/Award Number: CIFRE2020/0483; European Proteomics Infrastructure Consortium providing access, Grant/Award Number: 823839; H2020 Future and Emerging Technologies, Grant/Award Number: 829157; LABEX IBEID, Grant/Award Number: ANR-10-LABX-62-IBEID; Région Ile de France, Grant/Award Number: DIM 1HEALTH

Abstract

In antibody-based drug research, a complete characterization of antibody proteoforms covering both the amino acid sequence and all posttranslational modifications remains a major concern. The usual mass spectrometry-based approach to achieve this goal is bottom-up proteomics, which relies on the digestion of antibodies but does not allow the diversity of proteoforms to be assessed. Middle-down and top-down approaches have recently emerged as attractive alternatives but are not yet mastered and thus used in routine by many analytical chemistry laboratories. The work described here aims at providing guidelines to achieve the best sequence coverage for the fragmentation of intact light and heavy chains generated from a simple reduction of intact antibodies using Orbitrap mass spectrometry. Three parameters were found crucial to this aim: the use of an electron-based activation technique, the multiplex selection of precursor ions of different charge states, and the combination of replicates.

KEYWORDS

intact protein, LC-MS/MS, middle-down mass spectrometry, monoclonal antibodies, Orbitrap MS

1 | INTRODUCTION

In the last 20 years, antibodies (Abs) have become the new backbone of the pharmaceutical industry. The use of monoclonal Abs (mAbs) is now approved for over 67 targets, mostly involved in oncology but not only.^{1,2} Because mAbs are large molecules, their characterization requires the development of specific mass spectrometry approaches. Bottom-up proteomics (BUP), which is based on the digestion of Abs often with multiple enzymes, is the primary line of analysis, allowing both sequence information and posttranslational modification localization (PTMs) to be achieved in many cases.³⁻⁶ However, this

approach does not bring the required information on disulfide bond arrangements or simply on the different proteoforms present in a sample. These proteoforms can arise from N-terminal pyroglutamic acid formation, C-terminal clipping, glycosylation, glycation, methionine oxidation, deamidation, cysteinylolation, and so on.⁷ To fulfill this goal, alternative and complementary approaches have emerged, either on the intact mAbs (top-down proteomics or TDP) or on large parts of the Abs obtained either by S-S reduction or by specific enzymatic digestion⁸⁻¹¹ (middle-down proteomics or MDP). A recent interlaboratory study on the characterization of intact mAbs by TDP and MDP approaches including 20 partner laboratories and a wide

This is an open access article under the terms of the [Creative Commons Attribution-NonCommercial-NoDerivs](https://creativecommons.org/licenses/by-nc-nd/4.0/) License, which permits use and distribution in any medium, provided the original work is properly cited, the use is non-commercial and no modifications or adaptations are made.

© 2023 The Authors. *Journal of Mass Spectrometry* published by John Wiley & Sons Ltd.

diversity of instrumentation and methods highlighted the fact that there is no “one size fits all”.¹² One of the conclusions of this study is that compared with the techniques and instruments used for BUP, TDP and MDP show a clear potential for high-performance protein analysis but still require developments. One of the main advantages of TDP/MDP approaches is that the information achieved from MS/MS spectra is complementary to the BUP results and can also be obtained more rapidly. A very important piece of information is the molecular weight of intact mAbs. For extensive MS/MS characterization, the size and structure of these molecules (with both the intrachain and interchain S-S bonds) lead to incomplete sequence coverage that could be problematic to cover properly all complementarity determining regions (CDRs).¹³ Combining several activation techniques such as collision-based (collision-induced dissociation, CID; higher-energy collisional dissociation, HCD), electron-based (electron-transfer dissociation, ETD; electron-capture dissociation, ECD), combination of thereof (electron-transfer supplemented with collision or higher-energy collision dissociation, ETcId or ETHcD), or photon-based (ultra-violet photodissociation, UVPD) improves this issue but never leads to a 100% coverage.^{13–17} One strategy, developed by several groups, is to focus on the analysis of antibody subunits which can be either the light chain (Lc) and heavy chain (Hc) after disulfide reduction,^{18,19} Fab and Fc after specific enzymatic digestion^{13,20–22} or Lc, Fd and Fc/2 when combining both approaches.^{13,23–27} However, even with the smallest 25 kDa subunits a single MS/MS acquisition does not yield the expected sequence coverage and a combination of experiments are required. To date, no generic approach has been developed for this family of molecules, although such approach would be of high interest for biopharma companies. The objective of this paper is to provide a set of optimized instrumental parameters for the sequencing of Lc and Hc by Liquid Chromatography coupled to tandem Mass Spectrometry (LC–MS/MS) on an Orbitrap Tribrid mass spectrometer. Indeed, in contrast to BUP where the number of factors to optimize is very limited, in TDP/MDP, the quality of results strongly depends on a few parameters such as the precursor charge state,^{23,28–33} the activation technique and level (energy or time),^{6,13,18,19,23–27,29} and number of replicates. Compared with infusion experiments, the chromatographic separation brings additional constraints that will be discussed. Our results highlight the importance of electron-based fragmentation for improved sequence coverage. Initially restricted to a certain class of mass spectrometers (ECD in FT-ICR MS, ETD in ion traps), these methods are now more widely available on diverse platforms.

2 | MATERIAL AND METHODS

2.1 | Antibody sample preparation

Two commercially available mAbs (immunoglobulin type 1) were used for the fragmentation experiments: SigmaMab standard (SiLuLite, IgG1 lambda, CHO, Sigma) and NIST mAb standard (HzIgG1 kappa, NS0, NIST). These mAbs were chosen based on the strong sequence differences between their light chains (only 46% identity based on

sequence alignment). Each mAb was diluted to a final concentration of 0.32 µg/µL with guanidine hydrochloride (Sigma, 5 M final concentration) and reduced into Lc and Hc with DTT (Sigma, 100 mM final concentration) during 45 min at 45°C under 800 rpm. Samples were acidified with TFA (Sigma, 1% final concentration) before LC–MS/MS analysis.


2.2 | Liquid chromatography–mass spectrometry

LC–MS/MS analysis was performed using a Vanquish Horizon UHPLC system (Thermo Scientific, San Jose, CA) coupled to an Orbitrap Eclipse™ Tribrid mass spectrometer (Thermo Scientific, San Jose, CA) fitted with a H-ESI source. Ab1 µg of Lc and Hc was separated on a MAbPac™ RP column (2.1 mm i.d. × 50 mm, 4 µm particle size, 1500 Å pore size, Thermo Scientific, San Jose, CA) heated at 80°C with a gradient of 0.1% formic acid in acetonitrile (ACN) at 250 µL/min, from 25% to 41% in 2.8 min. For all experiments spray voltage was set to 3.8 kV, sheath gas settings was 35, auxiliary gas settings was 4, sweep gas settings was 1, vaporizer temperature was 150°C, ion transfer tube temperature was 350°C, RF value was 30%, and source fragmentation energy was 10 V. A first LC–MS experiment was acquired at 7500 resolving power (at m/z 200) with a scan range set to m/z 350–2500, five microscans (µscans) per MS scan, an automatic gain control (AGC) target value of 5×10^5 , and a maximum injection time of 50 ms. Fragmentation data were recorded using targeted LC–MS/MS experiments between 2.5 and 3.2 min for the NIST mAb Lc and 3.2 and 4.2 min for the NIST mAb Hc; between 2.6 and 3.25 min for the SigmaMab Lc and 3.2 and 4.0 min for the SigmaMab Hc. Four precursor charge states were chosen for each subunit across their respective charge state distribution, isolated by the quadrupole and subjected to individual or multiplexed fragmentation (Table S1). MS/MS scans were acquired at 120,000 resolving power (at m/z 200) with an isolation width of 1.6 m/z, five µscans, an AGC target value of 5×10^5 , and maximum injection time of 246 ms. HCD with normalized collision energies (NCE) of 15%, 20%, and 25%; ETHcD with 2, 5, and 10 ms of reaction time and a supplemental HCD activation with NCE of 5%, 10%, or 15%; and CID with collision energies (CE) of, respectively, 25%, 30%, 35% and 30%, 35%, 40% were used for the fragmentation of the Lc and Hc. For ETHcD experiments, the anionic fluoranthene reagent AGC target was set to 7×10^5 with a maximum injection time of 200 ms. All experiments were conducted using the Intact Protein mode with a pressure set to 1 mTorr in the ion-routing multipole (IRM).

2.3 | Data analysis

MS spectra were deconvolved with Genedata Expressionist® software using a time-resolved deconvolution and the Maximum Entropy (MaxEnt) algorithm. MS/MS spectra were averaged across the appropriate subunit elution windows and then deconvolved using the embedded Xtract algorithm in FreeStyle™ (v. 1.6.75.20) with a signal-

FIGURE 1 Chromatogram of NISTmAb Lc and Hc (top), MS spectra of Lc (bottom left), and Hc (bottom right). Precursor charge states selected for MS/MS analysis are depicted as red diamonds in MS spectra.

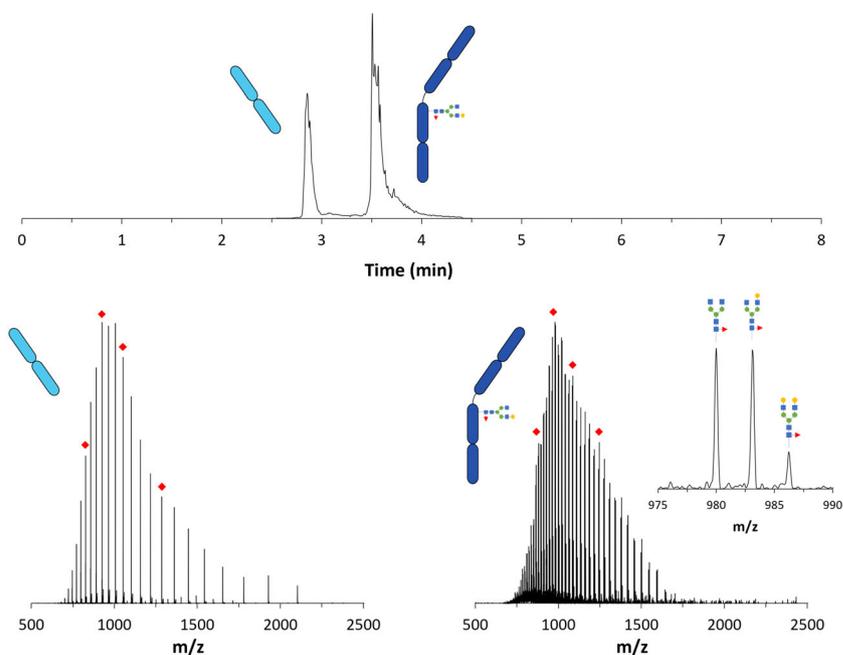

to-noise ratio threshold of 3, a fit factor of 80%, a remainder threshold of 25%, and maximum charge set to the precursor charge state. Lists of decharged and deisotoped monoisotopic masses were imported into ProSight Lite (v. 1.4)³⁴ and used for fragment assignments with a 5 ppm mass tolerance. Only b- and y-ions were considered for CID and HCD fragmentations, while b-, c-, y- and z-ions were searched for EThcD. TDFragMapper³⁵ was used for further visualization and comparison of fragmentation results. Lists of assigned fragments were exported from ProSight Lite and used with deconvolved data and the protein sequence as input for TDFragMapper. Venn diagrams and pairwise comparisons were generated using the Intervene Shiny app.³⁶ Combinations of fragmentation results from diverse MS/MS experiments were processed using an in-house developed R script.

3 | RESULT AND DISCUSSION

3.1 | Precursor charge state

Ionization of the reduced Lc and Hc from both NISTmAb and SigmaMab in denaturing conditions yields two charge state envelopes from 10+ to 30+ for the Lc and from 30+ to 65+ for the Hc (Figure 1 for NISTmAb and figure S1 for SigmaMab). From this envelope, a single charge state (ideally high charge with decent abundance) or several consecutive abundant charge states with a large isolation window are often selected for MS/MS fragmentation. Co-elution of multiple proteoforms in LC-MS can however hinder the use of a wide isolation window to avoid their simultaneous fragmentation and the production of chimeric spectra. In recent instruments, individual charge states can be either individually selected with a very narrow window or in a multiplexed manner to improve signal/noise ratio of the fragment ions. In top-down mass spectrometry, a clear link between the precursor

charge state and the quality of the MS/MS spectra has already been described and thus this parameter should be carefully optimized.^{23,28–33} To evaluate the role of this parameter when fragmenting Lc and Hc, four different charge states centered on the most abundant one in the envelope were selected either individually or simultaneously (multiplex). The lowest charge states (below 15+ for the Lc and below 40+ for the Hc) and highest ones were excluded from the study to avoid MS/MS spectra with very low intensity and competitive background noise. For the Hc, the three major glycoforms were co-selected for each charge state (Figure 1, right inset). As expected, very different raw and deconvolved MS/MS spectra were obtained for the HCD fragmentation of both NISTmAb Lc and Hc (Figure 2 and Figure S2). Similar results were obtained with other fragmentation techniques. As depicted on Figure 2, shifting from 17+ to 27+ (or multiplex) leads to very different fragmentation patterns, with a different distribution of fragment ions along both the mass and intensity axes. HCD activation of highly charged precursor ions generates a larger number of smaller fragments than the activation of lower charged precursors with the same activation energy. Fewer golden complementary pairs are also obtained when fragmenting highly charged precursors, probably due to the refragmentation of one of the ions. A golden complementary pair is a pair of fragment ions (a/x, b/y, or c/z) that has been formed by cleavage between the same pair of amino acids.³⁷ These golden complementary pairs are very useful for improving confidence in MS/MS sequencing. The smallest fragments (low mass region) can correspond to the cleavage of both the N- and C-terminal extremities and are therefore of interest to increase the sequence coverage (Figure S3). Although many of the observed fragments are common to all precursor charge states, unique ones can be retrieved for each individual precursor (both for the Lc and Hc subunits). Most of these unique fragments can be recovered in the multiplex experiments for the Lc (between 79% and 93%) but much less for

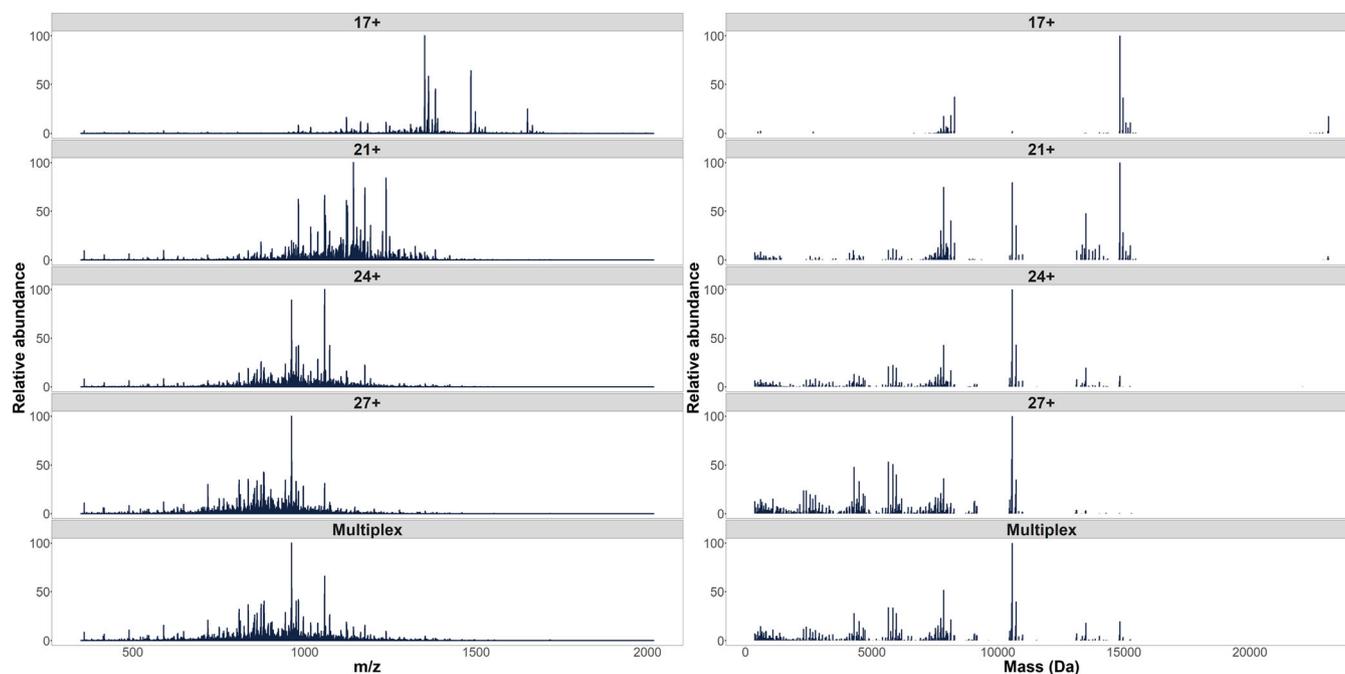


FIGURE 2 Higher-energy collisional dissociation (HCD) MS/MS spectra with NCE25% for each selected charge state of NISTmAb Lc before (left) and after deconvolution (right). From top to bottom: 17+, 21+, 24+, 27+, and multiplex.

TABLE 1 Comparison of residue cleavages between multiplex and the sum of individual precursors for the NISTmAb Lc and percentages of cleavages from individual fragmentations contained in multiplex experiments.

Fragmentation experiment	Cumulated residue cleavages (%)		Cleavages from individual precursors observed in multiplex (%)
	Multiplex	Individual precursors	
HCD 15%	30	33	86
HCD 20%	34	35	91
HCD 25%	34	37	87
ET2hcD5	65	77	80
ET2hcD15	63	77	79
ET5hcD5	74	79	90
ET5hcD15	75	74	92
ET10hcD10	73	66	93
ET10hcD15	68	58	92

Abbreviation: HCD, higher-energy collisional dissociation.

the Hc (between 30% and 90%) (Tables 1 and S2). Multiplexing several precursor charge states for the Lc analysis seems therefore to be a valuable strategy to combine most of the unique fragments of individual precursors into a single fragmentation experiment. This result obtained for a single replicate is also true for a triplicate: three multiplex experiments are almost equivalent to 12 individual ones (four precursors \times three replicates) for HCD but can also give better results for ETHcD (Table 1). For the Hc, the situation is slightly different and the sum of the 12 individual experiments always outperforms the multiplex ones for all tested HCD energies as well as for short ETHcD experiments (Table S2). Note however that we compare here the results of 12 experiments (individual precursors) to those obtained

only for three experiments (multiplex). Our results highlight that targeting different precursor charge states can provide complementary sequence information as different fragmentation pathways are activated. Moreover, using a multiplex selection is highly attractive as it reduces the number of experiments required to achieve the best sequence coverage.

3.2 | Fragmentation method

CID and HCD were compared with ETHcD (Figure 3). We decided not to include ETD because its efficiency as a standalone technique has

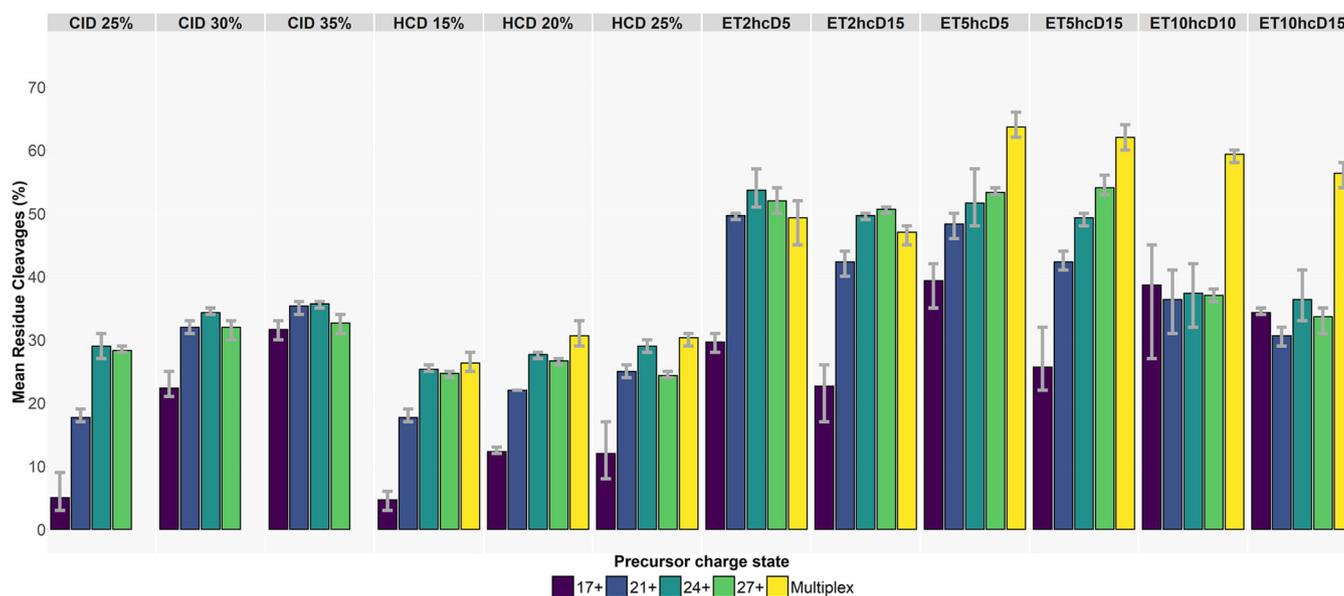


FIGURE 3 Mean residue cleavages (triplicates) for each selected precursor charge states of the NISTmAb Lc across the different fragmentation experiments.

already proved limited for the fragmentation of intact proteins, especially in an LC time scale. ETD is an attractive activation technique for fragmenting intact proteins but has to be combined with supplementary activation to yield satisfactory results and dissociate the fragments that are linked together by noncovalent interactions.^{30,38,39} The data obtained on the NISTmAb Lc indicate that both CID and HCD give similar results with a sequence coverage of 35% maximum. ETHcD gives the best results in all cases (up to 66% in a single experiment), independently of the precursor charge state studied here. Moreover, although multiplexing does not really bring a substantial improvement of CID/HCD, this is not the case for ETHcD especially with the largest activation times (5–10 ms). A similar conclusion can be drawn for the analysis of the NISTmAb Hc (Figure S4), with ETHcD leading to the best results despite a lower sequence coverage (maximum 19% for the multiplex in a single experiment) than what was achieved for the Lc. Similar results were obtained with the SigmaMab standard (Figures S5 and S6). We then compared the results obtained using each fragmentation method for the same precursor charge state (Figure 4 for NISTmAb Lc and S7 for NISTmAb Hc). The Venn diagrams confirm that ETHcD brings additional and unique cleavages compared with CID and HCD (although slower), with a charge state dependency. For the fragmentation of the NISTmAb Lc, the percentage of unique cleavages provided by ETHcD ranges from 37% to 46% for individual precursor ions and up to 67% for the multiplex. This is particularly the case for the higher charge states and less pronounced for the lowest ones for which CID leads to a significant number of unique cleavages. These unique cleavages represent 25% of all cleavages for the 17+ of the Lc but only 10% for the 40+ of the Hc, which is the lowest charge state fragmented. In contrast to CID, HCD on the Hc also brings unique cleavages (11%–18%) that can be explained by

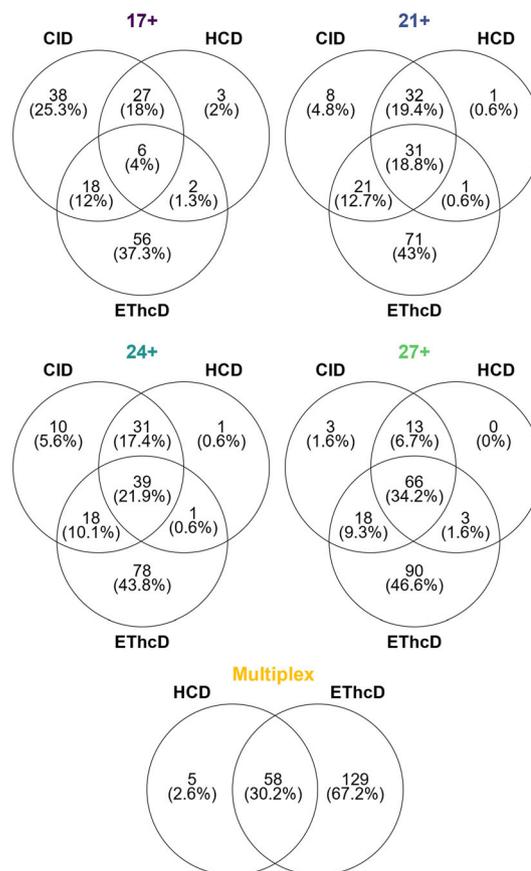


FIGURE 4 Venn diagrams of cleavages generated in collision-induced dissociation (CID), higher-energy collisional dissociation (HCD), and ETHcD for every targeted precursor charge state and multiplex for the NISTmAb Lc.

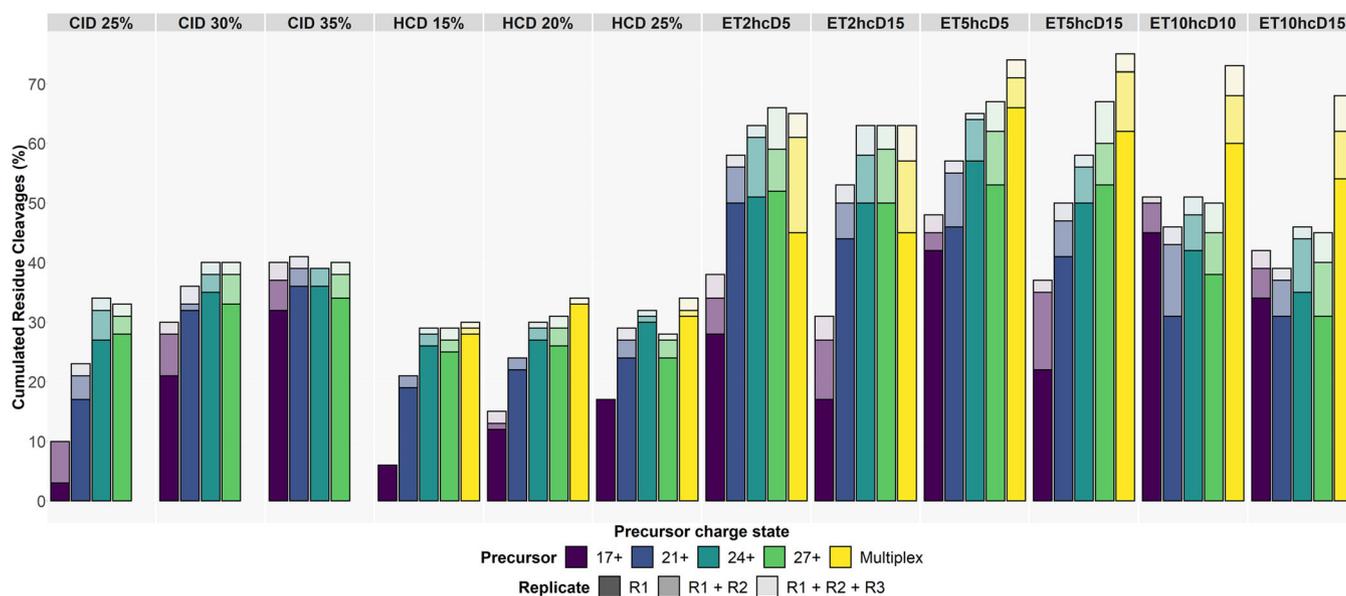


FIGURE 5 Summed residue cleavages over three replicates for each selected precursor charge states of the NISTmAb Lc across the different fragmentation experiments. R1, R2, and R3: Technical replicates numbered according to the order of acquisition (R1: first, R3: last).

the different energetic regime of both techniques. Concerning the general overlap of all three fragmentation techniques, the maximum achieved is 34% for the Lc and 8% for the Hc, which strongly supports the combination of fragmentation techniques for an improved sequence coverage. In conclusion, EThcD should always be favored since it leads to numerous unique fragments that cannot be obtained in HCD or CID.

3.3 | Technical replicates

Lc and Hc fragmentations take place on a chromatographic time scale, and therefore, the number of scans that can be accumulated is limited by the chromatographic peak width, which is 6 s in our conditions. Within this timeframe and with a resolution set to 120 k, a maximum of six scans of five microscans each can be acquired per LC run. Fragmentation of intact proteins with molecular masses above 25 kDa requires the accumulation of the largest number of scans to allow low intensity fragments to be extracted from the background noise. This scan accumulation remains one of the biggest challenges in LC-MS/MS experiments (compared with fragmentations acquired using infusion). One possible strategy to cope with this issue is to perform several replicates through repeated injections. Three replicates were therefore acquired for each fragmentation experiment in our study to assess the fragmentation repeatability and the expected added value of replicates that can account for signal variability⁴⁰ (example of the z_{106}^{9+} fragment ion of NISTmAb Lc is shown in Figure S8). Cumulated percentages of residue cleavages were computed and summarized for the NISTmAb Lc and Hc as shown in Figures 5 and S9, respectively. Our results indicate that increasing the number of replicates leads to a significant improvement for EThcD (up to 23% for the Lc and 11% for

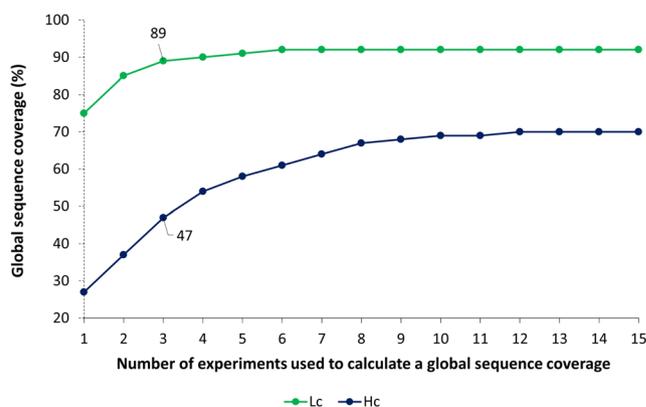


FIGURE 6 Maximum residue cleavage (%) obtained when combining diverse MS/MS experiments for the NISTmAb Lc and Hc.

the Hc with ET5hcD15) but much less for CID and HCD (maximum 9%/5% for Lc/Hc). The maximum sequence coverage obtained for three EThcD replicates is 75% for the Lc and 27% for the Hc, to be compared with 66% and 16%, respectively, for a single replicate. The improvement is therefore, as expected, more pronounced for the Hc. Similar observations could be made for the second mAb studied (Figures S10 and S11). The difference observed between EThcD and CID/HCD can probably be explained by the fact that EThcD opens more fragmentation channels than CID or HCD (with the formation of both *c/z* and *b/y*-ions and not only *b/y* ones), and thus, the fragment ions formed are of lower intensity. Cumulating replicates allows the diversity of EThcD cleavages to be better explored. It is also worth noting that the repeatability of the fragmentation is slightly influenced by the precursor charge state (and inherently by its relative intensity in the whole charge envelope). The conclusion here is that cumulating

at least two technical replicates, ideally three, improves the global sequence coverage of both Lc and Hc.

3.4 | Optimized combination of parameters

The objective here was to evaluate which combination of parameters, associated with the lowest number of experiments, can achieve the best sequence coverage for both Lc and Hc. In total, we performed 57 different experiments (in triplicate) for each subunit, varying the precursor charge state, the fragmentation type, and the activation level. The pairwise combination of all conditions would lead to 2^{57} possibilities (Figure S12), which is impossible to compute. We therefore decided to take the three best sets of parameters of each charge state and proceed with the comparison. These best sets are: two EThcD and one CID/HCD runs (detailed list in Table S3). The number of possibilities is now reduced to 2^{15} (32,768) which can be easily handled. The results are summarized in Figure 6 for NISTmAb (Figure S13 for SigmaMab). In line with the results described above, the combination of multiplex experiments always outperforms the others, whatever the types and level of fragmentation used, for both Lc and Hc. The best combination of three LC-MS/MS runs is composed by two EThcD experiments and one HCD. This shows that HCD still brings unique fragments that cannot be obtained by EThcD. These three best experiments lead to 89% sequence coverage for the NISTmAb Lc and 47% for the Hc. For the Lc, 89% is very close to the maximum of 92% achieved with six combined experiments (and more). These results are similar to the best ones obtained on the same Lc in the interlaboratory study¹² and better than the average sequence coverage achieved by all participants (50%). However, the situation is not as optimal for the Hc since the maximum is 70% for 12 and more combined experiments. As depicted in Figure 6, a plateau is rapidly reached for the Lc, although additional experiments always bring new and unique fragments for the Hc. This can be explained by the difference of size between the Lc and the Hc. It is well described that sequence coverage decreases with protein size.²² The sequence coverage we obtain for Hc of both mAbs (around 50% for three experiments) is much higher than the average observed (25%) in the interlaboratory study. Note however that the number of experiments summed in this study is unknown, which complicates a fair comparison. We then sought to make the same comparison but using only HCD, as the use of an unmodified benchtop Orbitrap Exploris would allow. The results are summarized in Figure S14 for both mAbs. The best combination of three LC-MS/MS runs with HCD activation corresponds to two multiplex and one experiment of the highest charge state and leads only to 42% for the Lc and 24% for the Hc. These results clearly show the added value of EThcD for improved sequence coverage.

4 | CONCLUSION

In this work, we extensively studied the combination of parameters leading to the best sequence coverage for the top-down analysis of

light and heavy chains formed from the reduction of intact Abs. It represents one of the most comprehensive studies performed on modern platforms offering experimental capabilities, such as multiplexing of charge states, previously not available. The parameters evaluated are the precursor charge state (individual and multiplex), the fragmentation type (HCD, CID, EThcD), and the fragmentation level (NCE for HCD, CE for CID, and reaction time and NCE for EThcD). In total, 171 LC-MS/MS runs were acquired (including triplicates) for each mAb: NISTmAb and SigmaMab. Similar behaviors were obtained for both Abs. Our results indicate that EThcD is mandatory to achieve the highest coverage of both Lc and Hc. HCD alone leads only to half cleavages, which clearly shows that Tribrid Orbitrap mass spectrometers remain to date the best platform for top-down MS/MS analysis. We also demonstrate that multiplexing leads to superior results than selecting individual precursor charge states. Multiplexing allows higher sequence coverage to be achieved in a lower number of experiments, which can be of interest in case of limited amount of sample. However, it is worth noting that, even with the best combination of parameters, a complete sequence coverage (100%) is never achieved, which leaves room for improvement. For example, the use of internal fragments could improve sequence coverage, in particular for CDRs, but only a few software tools have been developed to this aim⁴¹ and the confidence in their correct assignment can be challenging. However, the current trend toward an increased diversity of scaffolds in research pipelines of pharmaceutical companies will make this point even more critical in upcoming years.

AUTHOR CONTRIBUTIONS

The manuscript was written through contributions of all authors. All authors have given approval to the final version of the manuscript.

ACKNOWLEDGEMENTS

This work has received funding from the European Horizon 2020 Research and Innovation Program under the grant agreements no. 829157 (TopSpec) and 823839 (EPIC-XS). J. D acknowledges the Association Nationale de la Recherche et de la Technologie for funding his PhD under agreement CIFRE2020/0483. J.D. and J.C.R acknowledge the Region Ile de France (DIM 1HEALTH) and the LABEX IBEID (grant n° ANR-10-LABX-62-IBEID from the Programme d'Investissements d'Avenir) for funding the Eclipse Orbitrap.

DATA AVAILABILITY STATEMENT

Data are available via ProteomeXchange with identifier PXD038619.

ORCID

Jonathan Dhenin  <https://orcid.org/0000-0001-6615-2363>

REFERENCES

- Walsh G. Biopharmaceutical benchmarks 2018. *Nat Biotechnol.* 2018; 36(12):1136-1145. doi:10.1038/nbt.4305
- Kaplon H, Chenoweth A, Crescioli S, Reichert JM. Antibodies to watch in 2022. *MAbs.* 2022;14(1):2014296. doi:10.1080/19420862.2021.2014296

3. Millán-Martín S, Jakes C, Carillo S, et al. Inter-laboratory study of an optimised peptide mapping workflow using automated trypsin digestion for monitoring monoclonal antibody product quality attributes. *Anal Bioanal Chem.* 2020;412(25):6833-6848. doi:10.1007/s00216-020-02809-z
4. Liu T, Li Y, Xu J, et al. Mass spectrometry-based multi-attribute method for mutation analysis in the Early development of therapeutic proteins. *J Pharm Biomed Anal.* 2022;220:115018. doi:10.1016/j.jpba.2022.115018
5. Jakes C, Millán-Martín S, Carillo S, Scheffler K, Zaborowska I, Bones J. Tracking the behavior of monoclonal antibody product quality attributes using a multi-attribute method workflow. *J Am Soc Mass Spectrom.* 2021;32(8):1998-2012. doi:10.1021/jasms.0c00432
6. Dupré M, Duchateau M, Sternke-Hoffmann R, et al. De novo sequencing of antibody light chain proteoforms from patients with multiple myeloma. *Anal Chem.* 2021;93(30):10627-10634. doi:10.1021/acs.analchem.1c01955
7. Xu Y, Wang D, Mason B, et al. Structure, heterogeneity and developability assessment of therapeutic antibodies. *MAbs.* 2019;11(2):239-264. doi:10.1080/19420862.2018.1553476
8. Johansson BP, Shannon O, Björck L. IdeS: a bacterial proteolytic enzyme with therapeutic potential. *PLoS ONE.* 2008;3(2):e1692. doi:10.1371/journal.pone.0001692
9. von Pawel-Rammingen U, Johansson BP, Björck L. IdeS, a novel streptococcal cysteine proteinase with unique specificity for immunoglobulin G. *EMBO j.* 2002;21(7):1607-1615. doi:10.1093/emboj/21.7.1607
10. Moelleken J, Endesfelder M, Gassner C, et al. GingisKhan™ protease cleavage allows a high-throughput antibody to Fab conversion enabling direct functional assessment during lead identification of human monoclonal and bispecific IgG1 antibodies. *MAbs.* 2017;9(7):1076-1087. doi:10.1080/19420862.2017.1364325
11. Sjögren J, Olsson F, Beck A. Rapid and improved characterization of therapeutic antibodies and antibody related products using IdeS digestion and subunit analysis. *Analyst.* 2016;141(11):3114-3125. doi:10.1039/C6AN00071A
12. Srzentić K, Fornelli L, Tsybin YO, et al. Interlaboratory study for characterizing monoclonal antibodies by top-down and middle-down mass spectrometry. *J Am Soc Mass Spectrom.* 2020;31(9):1783-1802. doi:10.1021/jasms.0c00036
13. Fornelli L, Srzentić K, Huguet R, et al. Accurate sequence analysis of a monoclonal antibody by top-down and middle-down Orbitrap mass spectrometry applying multiple ion activation techniques. *Anal Chem.* 2018;90(14):8421-8429. doi:10.1021/acs.analchem.8b00984
14. Mao Y, Valeja SG, Rouse JC, Hendrickson CL, Marshall AG. Top-down structural analysis of an intact monoclonal antibody by electron capture dissociation-Fourier transform ion cyclotron resonance-mass spectrometry. *Anal Chem.* 2013;85(9):4239-4246. doi:10.1021/ac303525n
15. Lodge JM, Schauer KL, Brademan DR, et al. Top-down characterization of an intact monoclonal antibody using activated ion electron transfer dissociation. *Anal Chem.* 2020;92(15):10246-10251. doi:10.1021/acs.analchem.0c00705
16. Tsybin YO, Fornelli L, Stoermer C, et al. Structural analysis of intact monoclonal antibodies by electron transfer dissociation mass spectrometry. *Anal Chem.* 2011;83(23):8919-8927. doi:10.1021/ac201293m
17. Fornelli L, Damoc E, Thomas PM, et al. Analysis of intact monoclonal antibody IgG1 by electron transfer dissociation Orbitrap FTMS. *Mol Cell Proteomics.* 2012;11(12):1758-1767. doi:10.1074/mcp.M112.019620
18. Melani RD, Srzentić K, Gerbası VR, et al. Direct measurement of light and heavy antibody chains using ion mobility and middle-down mass spectrometry. *MAbs.* 2019;11(8):1351-1357. doi:10.1080/19420862.2019.1668226
19. He L, Anderson LC, Barnidge DR, Murray DL, Hendrickson CL, Marshall AG. Analysis of monoclonal antibodies in human serum as a model for clinical monoclonal gammopathy by use of 21 tesla FT-ICR top-down and middle-down MS/MS. *J Am Soc Mass Spectrom.* 2017;28(5):827-838. doi:10.1007/s13361-017-1602-6
20. Greisch J-F, den Boer MA, Beurskens F, et al. Generating informative sequence tags from antigen-binding regions of heavily glycosylated IgA1 antibodies by native top-down electron capture dissociation. *J Am Soc Mass Spectrom.* 2021;32(6):1326-1335. doi:10.1021/jasms.0c00461
21. Bondt A, Hoek M, Tamara S, et al. Human plasma IgG1 repertoires are simple, unique, and dynamic. *Cell Syst.* 2021;12(12):1131-1143. e5. doi:10.1016/j.cels.2021.08.008
22. Srzentić K, Nagornov KO, Fornelli L, et al. Multiplexed middle-down mass spectrometry as a method for revealing light and heavy chain connectivity in a monoclonal antibody. *Anal Chem.* 2018;90(21):12527-12535. doi:10.1021/acs.analchem.8b02398
23. Shaw JB, Malhan N, Vasil'ev YV, et al. Sequencing grade tandem mass spectrometry for top-down proteomics using hybrid electron capture dissociation methods in a benchtop Orbitrap mass spectrometer. *Anal Chem.* 2018;90(18):10819-10827. doi:10.1021/acs.analchem.8b01901
24. Fornelli L, Ayoub D, Aizikov K, Beck A, Tsybin YO. Middle-down analysis of monoclonal antibodies with electron transfer dissociation Orbitrap Fourier transform mass spectrometry. *Anal Chem.* 2014;86(6):3005-3012. doi:10.1021/ac4036857
25. Watts E, Williams JD, Miesbauer LJ, Bruncko M, Brodbelt JS. Comprehensive middle-down mass spectrometry characterization of an antibody-drug conjugate by combined ion activation methods. *Anal Chem.* 2020;92(14):9790-9798. doi:10.1021/acs.analchem.0c01232
26. Hernandez-Alba O, Houel S, Hessmann S, et al. A case study to identify the drug conjugation site of a site-specific antibody-drug-conjugate using middle-down mass spectrometry. *J Am Soc Mass Spectrom.* 2019;30(11):2419-2429. doi:10.1007/s13361-019-02296-2
27. Cotham VC, Brodbelt JS. Characterization of therapeutic monoclonal antibodies at the subunit-level using middle-down 193 nm ultraviolet photodissociation. *Anal Chem.* 2016;88(7):4004-4013. doi:10.1021/acs.analchem.6b00302
28. Brunner AM, Lössl P, Liu F, et al. Benchmarking multiple fragmentation methods on an orbitrap fusion for top-down phospho-proteome characterization. *Anal Chem.* 2015;87(8):4152-4158. doi:10.1021/acs.analchem.5b00162
29. Shliaha PV, Gibb S, Gorshkov V, et al. Maximizing sequence coverage in top-down proteomics by automated multimodal gas-phase protein fragmentation. *Anal Chem.* 2018;90(21):12519-12526. doi:10.1021/acs.analchem.8b02344
30. Riley NM, Westphall MS, Coon JJ. Sequencing larger intact proteins (30-70 KDa) with activated ion electron transfer dissociation. *J Am Soc Mass Spectrom.* 2018;29(1):140-149. doi:10.1007/s13361-017-1808-7
31. Riley NM, Mullen C, Weisbrod CR, et al. Enhanced dissociation of intact proteins with high capacity electron transfer dissociation. *J Am Soc Mass Spectrom.* 2016;27(3):520-531. doi:10.1007/s13361-015-1306-8
32. Durbin KR, Fellers RT, Ntai I, Kelleher NL, Compton PD. Autopilot: an online data acquisition control system for the enhanced high-throughput characterization of intact proteins. *Anal Chem.* 2014;86(3):1485-1492. doi:10.1021/ac402904h
33. Shaw JB, Li W, Holden DD, et al. Complete protein characterization using top-down mass spectrometry and ultraviolet photodissociation. *J Am Chem Soc.* 2013;135(34):12646-12651. doi:10.1021/ja4029654
34. Fellers RT, Greer JB, Early BP, et al. ProSight lite: graphical software to analyze top-down mass spectrometry data. *Proteomics.* 2015;15(7):1235-1238. doi:10.1002/pmic.201570050

35. Dhenin J, Borges Lima D, Dupré M, Chamot-Rooke J. TDFragMapper: a visualization tool for evaluating experimental parameters in top-down proteomics. *Bioinformatics*. 2021;38(4):1136-1138. doi:10.1093/bioinformatics/btab784
36. Khan A, Mathelier A. Intervene: a tool for intersection and visualization of multiple gene or genomic region sets. *BMC Bioinform*. 2017;18(1):287. doi:10.1186/s12859-017-1708-7
37. Horn DM, Zubarev RA, McLafferty FW. Automated de novo sequencing of proteins by tandem high-resolution mass spectrometry. *Proc Natl Acad Sci*. 2000;97(19):10313-10317. doi:10.1073/pnas.97.19.10313
38. Riley NM, Sikora JW, Seckler HS, et al. The value of activated ion electron transfer dissociation for high-throughput top-down characterization of intact proteins. *Anal Chem*. 2018;90(14):8553-8560. doi:10.1021/acs.analchem.8b01638
39. Riley NM, Westphall MS, Coon JJ. Activated ion-electron transfer dissociation enables comprehensive top-down protein fragmentation. *J Proteome Res*. 2017;16(7):2653-2659. doi:10.1021/acs.jproteome.7b00249
40. Wiener MC, Sachs JR, Deyanova EG, Yates NA. Differential mass spectrometry: a label-free LC-MS method for finding significant differences in complex peptide and protein mixtures. *Anal Chem*. 2004;76(20):6085-6096. doi:10.1021/ac0493875
41. Lantz C, Zenaidee MA, Wei B, Hemminger Z, Ogorzalek Loo RR, Loo JA. ClipsMS: an algorithm for analyzing internal fragments resulting from top-down mass spectrometry. *J Proteome Res*. 2021;20(4):1928-1935. doi:10.1021/acs.jproteome.0c00952

SUPPORTING INFORMATION

Additional supporting information can be found online in the Supporting Information section at the end of this article.

How to cite this article: Dhenin J, Dupré M, Druart K, Krick A, Mauriac C, Chamot-Rooke J. A multiparameter optimization in middle-down analysis of monoclonal antibodies by LC-MS/MS. *J Mass Spectrom*. 2023;58(3):e4909. doi:10.1002/jms.4909